

# Purification Strategies for Flavones

and Related Compounds  
with CombiFlash Systems

## Abstract

CombiFlash<sup>®</sup> purification systems, in conjunction with RediSep<sup>®</sup> columns, are useful for purifying flavonoid compounds. Examples are provided from the scientific literature. Flavonoid compounds are purified on silica using hexane/ethyl acetate or dichloromethane/methanol gradients. Reverse phase C18 is a useful method to purify this class of compounds.

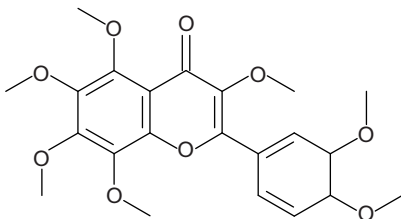
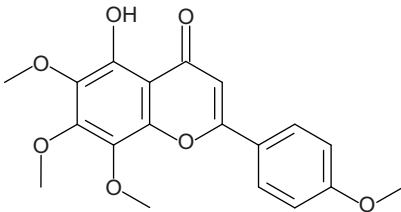
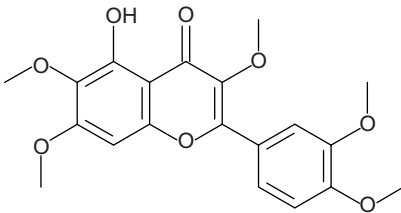
## Introduction

Flavonoids are bioactive compounds ubiquitous to green plants. These compounds are often colored and show anti-inflammatory, anti-carcinogenic, and anti-atherogenic properties<sup>1</sup>. The unique medicinal properties of these compounds make them attractive candidates for purification by natural products chemists and synthesis by medicinal chemists.

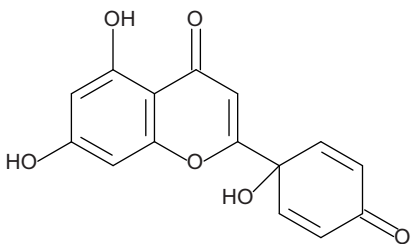
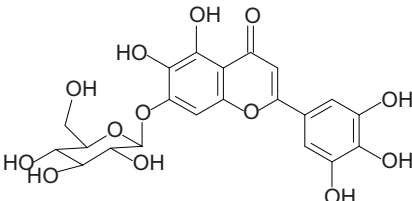
## Purification Strategies

Table 1 shows that flavonoid compounds are moderately polar. When the phenol groups are methylated (compounds **1** to **3**), the methoxyflavones are purified on silica with a hexane-ethyl acetate gradient. As the number of phenol groups are increased, the increased polarity of the flavones require more polar solvents such as methanol (compound **4**). A flavone (compound **5**) with a glycoside showing further increases in polarity was purified on C18 rather than silica. Typical absorption (detection) wavelengths for these compounds range between 250–285 nm and also 320–380 nm<sup>2</sup>. The purifications listed showed no use of organic modifiers. When running reverse phase, the use of 0.1% trifluoroacetic acid, formic acid, or acetic acid may be useful to sharpen peaks. When purifying the compounds on silica gel, using an ethyl acetate/methanol gradient may prove useful as ethyl acetate is in a different selectivity group than the chlorinated solvents<sup>3</sup>.

**Table 1: Purification of non-polar flavones**

ID and Reference	Structure	Name	Solvent System
1 <sup>4</sup>		3',4',3,5,6,7,8-Heptamethoxyflavone	Silica column: hexane: ethyl acetate 33 to 50% ethyl acetate gradient  C18 column: water: methanol 20 to 80% methanol gradient
2 <sup>2</sup>		5-hydroxy-6,7,8,4'-tetramethoxyflavone	Silica column: hexane/ethyl acetate gradient
3 <sup>1</sup>		5-hydroxy-3,6,7,3',4'-pentamethoxyflavone	Silica column: hexane/ethyl acetate gradient

**Table 1: Purification of non-polar flavones**

ID and Reference	Structure	Name	Solvent System
4 <sup>5</sup>		Protoapigenone	33% MeOH: 67% CHCl <sub>3</sub> isocratic
5 <sup>6</sup>		5,6,7,3',4',5'-hexahydroxyflavone-7-O-β-glucopyranoside	C18 column: water: methanol gradient

## Conclusion

Flavonoid compounds are easily purified with CombiFlash systems with RediSep columns. The compounds exhibit modest polarity and can be purified with hexane/ethyl acetate or dichloromethane/methanol gradients. Reverse phase chromatography can also be employed to purify this compound class.

## References

- Li, S.; Pan, M-H; Lai, C-S; Lo, C-Y; Dushenkova, S.; Ho, C-T. Isolation and syntheses of polymethoxyflavones and hydroxylated polymethoxyflavones as inhibitors of HL-60 cell lines. *Bioorg. Med. Chem.* **2007**, *15*, 3381–3389.
- Li, S.; Lo, C-Y; Ho, C-T Hydroxylated Polymethoxyflavones and Methylated Flavonoids in Sweet Orange (*Citrus sinensis*) Peel. *J. Agric. Food Chem.* **2006**, *54*, 4716–4185.
- Effective Organic Compound Purification*. 4<sup>th</sup> Ed., Teledyne Isco, Lincoln NE, 2010, 105.
- Manthey, J.A.; Bendele, P. Anti-inflammatory Activity of an Orange Peel Polymethoxylated Flavone, 3',4',3,5,6,7,8-Heptamethoxyflavone, in the Rat Carrageenan/Paw Edema and Mouse Lipopolysaccharide-Challenge Assays. *J. Agric. Food Chem.* **2008**, *56* (20), 9399-9403.
- Lin, A-S; Nakagawa-Goto, K.; Chang, F-R; Yu, D.; Morris-Natschke, S.L.; Wu, C-C; Chen, S-L; Yang-Chang Wu, Y-C; Lee, K-H. First Total Synthesis of Protoapigenone and Its Analogues as Potent Cytotoxic Agents. *J. Med. Chem.* **2007**, *50*, 3921-3927.
- Meng, Y.; Krzysiak, A.J.; Durako, M.J.; Kunzelman, J.I.; L.C. Wright, J.L.C. Flavones and flavone glycosides from *Halophila johnsonii*. *Phytochemistry*, **2008**, *69* (14), 2603-2608.

Last modified November 14, 2011

### Teledyne Isco, Inc.

P.O. Box 82531, Lincoln, Nebraska, 68501 USA  
Toll-free: (800) 775-2965 • Phone: (402) 464-0231 • Fax: (402) 465-3001  
E-mail: IscoService@teledyne.com

Teledyne Isco is continually improving its products and reserves the right to change product specifications, replacement parts, schematics, and instructions without notice.

